

# ACCURUN® 52 Multi-Marker Hepatitis Positive Control 2

## NAME AND INTENDED USE

ACCURUN® controls are intended for use as positive qualitative controls to monitor laboratory testing precision and detect errors in laboratory testing procedures. ACCURUN 52 Multi-Marker Hepatitis Positive Control 2 has been formulated for use with *in vitro* diagnostic test kits for the qualitative determination of antibodies to Hepatitis B e Antigen (anti-HBe), Hepatitis B Surface Antigen (anti-HBs), and Hepatitis A Virus (anti-HAV). This product is not intended for use in testing blood or plasma donors. ACCURUN controls do not have quantitative assigned values. *For professional laboratory use only.*

## SUMMARY

Frequent testing of independent quality control samples provides the analyst with a means of monitoring the performance of laboratory assays. Routine use of controls enables laboratories to monitor day-to-day test variation, lot-to-lot performance of test kits, and operator variation, and can assist in identifying increases in random or systematic error. A well designed quality control program can provide added confidence in the reliability of results obtained for unknown specimens. The use of low-reactive samples as independent controls may provide valuable information concerning laboratory proficiency and kit lot variation that may affect assay sensitivity<sup>1</sup>.

## PRINCIPLES OF THE PROCEDURE

ACCURUN 52 Multi-Marker Hepatitis Positive Control 2 has been designed for use with *in vitro* assay procedures for purposes of monitoring test performance. ACCURUN 52 Multi-Marker Hepatitis Positive Control 2 is manufactured from human serum or plasma reactive for anti-HBe, anti-HBs and anti-HAV, and nonreactive for HBsAg and antibodies to HIV 1 and 2, HTLV, and HCV. ACCURUN controls do not have assigned values. Examples of assay(s) with which this control may be compatible are listed in Table 1. Specific levels of reactivity will vary among different manufacturers' assays, different procedures, different lot numbers, and different laboratories.

## REAGENTS

Item No. 2000-0043

6 vials, 3.5 ml per vial

This control contains stabilizers (EDTA, buffering agents) and 0.1% ProClin® (5-chloro-2-methyl-4-isothiazolin-3-one & 2-methyl-4-isothiazolin-3-one) as preservative.

## WARNINGS AND PRECAUTIONS

### For In Vitro Diagnostic Use

**CAUTION:** Handle ACCURUN controls and all human blood products as though capable of transmitting infectious agents. ACCURUN 52 Multi-Marker Hepatitis Positive Control 2 is manufactured from human serum or plasma reactive for anti-HBe, anti-HBs, and anti HAV. The serum or plasma used in manufacturing is nonreactive for HBsAg and antibodies to HIV 1 and 2, HTLV, and HCV with current FDA licensed tests.

### Safety Precautions

Use the Centers for Disease Control (CDC) recommended universal precautions for handling ACCURUN controls and human blood<sup>2</sup>. Do not pipette by mouth; do not eat or drink in areas where specimens are being handled. Clean any spillage by immediately wiping up with 0.5% sodium hypochlorite solution. Dispose of all specimens, controls and materials used in testing as though they contain infectious agents. Additional safety information can be found in the product Safety Data Sheet (SDS) found on the company website.

### Handling Precautions

Do not use ACCURUN controls beyond the expiration date. Avoid microbial contamination of the controls when opening and closing the vials.

## STORAGE INSTRUCTIONS

Store ACCURUN controls refrigerated at 2-8°C. Once opened, ACCURUN controls should be stored at 2-8°C when not in use and discarded after 60 days. After opening, record the date opened and the expiration date on the vial. Multiple freeze-thaw cycles are not recommended and may have variable adverse effects upon test results. To prevent leakage, store vials upright.

## INDICATIONS OF REAGENT INSTABILITY OR DETERIORATION

Alterations in physical appearance may indicate instability or deterioration of ACCURUN controls. Solutions that are visibly turbid should be discarded.

## PROCEDURE

### Materials Provided

ACCURUN 52 Multi-Marker Hepatitis Positive Control 2 is manufactured from human serum or plasma, including materials reactive for anti-HBe, anti-HBs and anti-HAV and nonreactive for HBsAg and antibodies to HIV 1 and 2, HTLV and HCV. See REAGENTS for a list of package sizes.

### Materials Required but not Provided

Refer to instructions supplied by manufacturers of the test kits to be used.

### Instructions for Use

Allow the controls to reach room temperature prior to use, then return controls to refrigerated storage immediately after use. Mix the contents of the vials by gently swirling. ACCURUN controls should be included in a test run using exactly the same procedure provided by the manufacturer for unknown specimens. ACCURUN controls must NOT be substituted for the positive and negative control reagents provided with licensed test kits.

## Quality Control

Since ACCURUN controls do not have assigned values, it is recommended that each laboratory validate the use of each lot of ACCURUN control with each specific assay system prior to its routine use in the laboratory.

## INTERPRETATION OF RESULTS

Levels of reactivity of ACCURUN 52 Multi-Marker Hepatitis Positive Control 2 may vary with different manufacturers' tests and different test kit lots. Each laboratory must establish a range of acceptable values for ACCURUN Controls. When results for ACCURUN 52 Multi-Marker Hepatitis Positive Control 2 are outside the established acceptable range of values, it may be an indication of unsatisfactory test performance. Possible sources of error are: deterioration of test kit reagents, operator error, faulty performance of equipment, or contamination of reagents.

## LIMITATIONS OF THE PROCEDURE

ACCURUN CONTROLS MUST NOT BE SUBSTITUTED FOR THE POSITIVE AND NEGATIVE CONTROL REAGENTS PROVIDED WITH MANUFACTURED TEST KITS.

TEST PROCEDURES and INTERPRETATION OF RESULTS provided by manufacturers of test kits must be followed closely. Deviations from procedures recommended by test kit manufacturers may produce unreliable results. ACCURUN controls are qualitative, not automated, and are provided for quality assurance purposes and must not be used for calibration or as a primary reference preparation in any test procedure. Adverse shipping and/or storage conditions or use of outdated controls may produce erroneous results.

## EXPECTED RESULTS

**ACCURUN Controls DO NOT HAVE AN ASSIGNED VALUES.** This control has been formulated to produce positive reactivity for the analytes listed in Table 1. Specific levels of reactivity will vary among different manufacturers' assays, different procedures, different reagent lot numbers, and different laboratories. Each laboratory should establish its own range of acceptable values for each analyte. For example, the acceptable range might include all values within 2 standard deviations of the mean of 20 data points obtained in 20 runs over a period of 30 days<sup>3</sup>.

## SPECIFIC PERFORMANCE CHARACTERISTICS

ACCURUN controls have been designed for use with *in vitro* assay procedures for purposes of monitoring assay performance. ACCURUN 52 Multi-Marker Hepatitis Positive Control 2 is manufactured from human serum or plasma including materials reactive for anti-HBe, anti-HBs and anti-HAV and nonreactive for antibodies to HIV 1 and 2, HTLV and HCV. ACCURUN controls do not have assigned values. Specific levels of reactivity will vary among different manufacturers' assays, different procedures, different reagent lot numbers, and different laboratories. Procedures for implementing a quality assurance program and monitoring test performance on a routine basis must be established by each individual laboratory. Quality control materials should be used in accordance with local, state, and federal regulations and accreditation requirements.

## REFERENCES

- Green IV GA, Carey RN, Westgard JO, Carten T, Shablesky LA, Achord D, Page E, and Le AV. Quality control for qualitative assays: quantitative QC procedure designed to assure analytical quality required for an ELISA for hepatitis B surface antigen. Clin. Chem. 43:9 1618-1621, 1997.
- Siegel JD, Rhinehart E, Jackson M, Chiarello L, and the Healthcare Infection Control Practices Advisory Committee, 2007 Guideline for Isolation Precautions: Preventing Transmission of Infectious Agents in Healthcare Settings.
- Statistical Quality Control for Quantitative Measurements: Principles and Definitions; Approved Guideline— Fourth Edition. CLSI document C24, 2016.

Table 1. This product is tested at release using the following assays.

Marker	Manufacturer/Product Name
Anti-HBe	DiaSorin LIAISON® Anti-HBe
Anti-HBs	DiaSorin LIAISON® Anti-HBs
Anti-HAV	Bio-Rad MONOLISA anti-HAV EIA


For assistance, contact LGC Clinical Diagnostics Technical Support at +1 508.244.6400.

Any serious incident that has occurred in relation to the device shall be reported to LGC Clinical Diagnostics Technical Support and, if in use in the EU, the competent authority of the Member State in which the incident occurred.

Date	Description of Change
May 2025	Update for IVDR

# ACCURUN® 52

Multi-Marker Hepatitis  
Positive Control 2

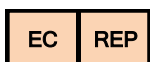

MEDIMARK® Europe  
11, rue Émile Zola BP 2332  
38033 Grenoble Cedex 2 – France  
+ 33 (0) 4 76 86 43 22  
info@medimark-europe.com



LGC Clinical Diagnostics, Inc. | 37 Birch Street, Milford, MA 01757 USA  
Phone: +1 508.244.6400 | CDx-Info@LGCGroup.com

10127en-GB-13 May 2025

## Legend of all symbols used in LGC product labeling



Upper limit of temperature



Temperature limitation



Authorized Representative in  
the European Community



Biological risks



Use By



*In Vitro* Diagnostic Medical Device



Negative control



Catalogue number



Consult instructions for use



Positive control



Batch code



Manufacturer



Control



Highly Flammable



Toxic by inhalation, in contact  
with skin and if swallowed



Health Hazard



Importer