


ACCURUN®

HIV RNA
Tri-Level Controls Kit

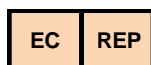

SeraCare Life Sciences, Inc. | 37 Birch Street, Milford, MA 01757 USA
Phone: +1 508.244.6400 | info@seracare.com

MEDIMARK® Europe
11, rue Émile Zola BP 2332
38033 Grenoble Cedex 2 – France
+ 33 (0) 4 76 86 43 22
info@medimark-europe.com

12967GB-02 May 2017

Explanation of symbols used in SeraCare product labeling



Upper limit of temperature



Temperature limitation



Authorized Representative in
the European Community



Biological risks



Use By



In Vitro Diagnostic Medical Device



Negative control



Catalogue number



Consult instructions for use



Positive control



Batch code



Manufacturer



Control



Highly Flammable



Toxic by inhalation, in contact
with skin and if swallowed



Health Hazard

ACCURUN® HIV RNA Tri-Level Controls Kit

NAME AND INTENDED USE

ACCURUN products are intended to estimate laboratory testing precision and can be used to detect errors in laboratory testing procedures. ACCURUN® HIV RNA Tri-Level Controls Kit is formulated for use with *in vitro* diagnostic amplified test methods that detect and quantitate Human Immunodeficiency Virus (HIV) RNA. *For In Vitro Diagnostic Use.*

SUMMARY

Frequent testing of independent quality control samples provides the analyst with a means of monitoring the performance of laboratory assays. Routine use of controls enables laboratories to detect immediate analytical errors and monitor long term performance of test kits, and can assist in identifying increases in random or systematic error. A well-designed quality control program can provide added confidence in the reliability of results obtained for unknown specimens. The use of independent controls may provide valuable information concerning laboratory proficiency and kit lot variation that may affect assay sensitivity¹.

PRINCIPLES OF THE PROCEDURE

ACCURUN HIV RNA Tri-Level Controls Kit is designed for use with *in vitro* diagnostic assay procedures for the purpose of monitoring test performance. This product is a set of external run controls for use with assays that detect and quantitate HIV. ACCURUN HIV RNA Tri-Level Controls kit is manufactured from defibrinated human plasma and HIV RNA and has been formulated for positive and negative reactivity with quantitative HIV RNA assays. The kit contains three controls: a high positive, a low positive, and a negative control.

REAGENTS

Item No. 2020-0142	15 vials/kit
	5 vials of each control:
	HIV RNA Tri-Level High, 1.2 ml per vial
	HIV RNA Tri-Level Low, 1.2 ml per vial
	HIV RNA Tri-Level Neg, 1.2 ml per vial

ACCURUN HIV RNA Tri-Level Controls Kit contains stabilizers and 0.09% sodium azide as a preservative.

WARNINGS AND PRECAUTIONS

For In Vitro Diagnostic Use

CAUTION: Handle ACCURUN controls and all human derived products as though capable of transmitting infectious agents. ACCURUN HIV RNA Tri-Level Controls Kit is manufactured using defibrinated human plasma and HIV RNA positive plasma.

Safety Precautions

Use Center for Disease Control and Prevention (CDC) recommended universal precautions for handling ACCURUN and human specimens². Do not pipette by mouth; do not smoke, eat or drink in areas where specimens are being handled. Clean any spillage by immediately wiping up with 0.5% sodium hypochlorite solution. Dispose of all controls, specimens and materials used in testing as though they contain infectious agents³.

Handling Precautions

Do not use ACCURUN controls beyond the expiration date. Avoid contamination of the controls when opening and closing the vials.

STORAGE INSTRUCTIONS

Store ACCURUN HIV RNA Tri-Level Controls Kit at -70 °C or colder until use. Once opened, an individual vial of ACCURUN Tri-Level Control should not be reused. Store vials upright to prevent leakage.

INDICATIONS OF REAGENT INSTABILITY OR DETERIORATION

Alterations in physical appearance may indicate instability or deterioration of ACCURUN Controls. Solutions that are visibly turbid should be discarded.

PROCEDURE

Materials Provided

ACCURUN HIV RNA Tri-Level Controls kit is manufactured from defibrinated human plasma and HIV RNA. It also contains human proteins and preservatives.

Materials Required but not Provided

Refer to instructions supplied by manufacturers of the test kits to be used.

Instructions for Use

Allow the controls to come to room temperature and mix the contents of the vials by gentle inversion before use. Once thawed, use immediately. ACCURUN controls should be included in a test run using the same procedure provided by the manufacturer for unknown specimens. ACCURUN controls must NOT be substituted for the positive and negative control reagents provided with the manufactured test kits.

Quality Control

Since ACCURUN HIV RNA Tri-Level Controls do not have assigned values, it is recommended that each laboratory validate the use of each lot of ACCURUN with each specific assay system prior to its routine use in the laboratory.

INTERPRETATION OF RESULTS

Levels of reactivity of ACCURUN HIV RNA Tri-Level Controls may vary with different manufacturers' tests and different test kit lots. Since the controls do not have assigned values, the laboratory must establish a range for each lot of each control. When results for ACCURUN HIV RNA Tri-Level Controls are outside of the established acceptance range for any of the three controls, it may indicate unsatisfactory test performance. Possible sources of error include: deterioration of test kit reagents, operator error, faulty performance of equipment, or contamination of reagents.

LIMITATIONS OF THE PROCEDURE

ACCURUN CONTROLS MUST NOT BE SUBSTITUTED FOR ANY CONTROL REAGENTS PROVIDED WITH MANUFACTURED TEST KITS.

TEST PROCEDURES and INTERPRETATION OF RESULTS provided by manufacturers of test kits must be followed closely. Deviations from procedures recommended by test kit manufacturers may produce unreliable results. ACCURUN controls are not calibrators and should not be used for assay calibration. Performance characteristics for ACCURUN HIV RNA Tri-Level Controls Kit have been established only for amplified nucleic acid tests that detect HIV RNA. Adverse shipping and storage conditions or use of outdated product may produce erroneous results.

EXPECTED RESULTS

ACCURUN HIV RNA Tri-Level Controls Kit DOES NOT HAVE ASSIGNED VALUES. Specific levels of reactivity will vary among different manufacturers' assays, different procedures, different lot numbers, and different laboratories. ACCURUN HIV RNA Tri-Level Controls Kit contains a high positive control, a low positive control, and a negative control. Each control is formulated to produce different results. Each laboratory should establish its own range of acceptable values for each individual control. Refer to Table 1 for typical results for each control in the ACCURUN HIV RNA Tri-Level Controls Kit.

SPECIFIC PERFORMANCE CHARACTERISTICS

ACCURUN controls have been designed for use with *in vitro* diagnostic assay procedures for the purposes of monitoring assay performance. ACCURUN controls do not have assigned values. Procedures for implementing a quality assurance program and monitoring test performance on a routine basis must be established by each individual laboratory.

REFERENCES

- Green IV GA, Carey RN, Westgard JO, Carten T, Shablesky LA, Achord D, Page E, and Le AV. Quality control for qualitative assays: quantitative QC procedures designed to assure analytical quality required for an ELISA for hepatitis B surface antigen. Clin. Chem. 43:9 1618-1621, 1997.
- Siegel JD, Rhinehart E, Jackson M, Chiarello L, and the Healthcare Infection Control Practices Advisory Committee, 2007 Guideline for Isolation Precautions: Preventing Transmission of Infectious Agents in Healthcare Settings.
- Treatment Standards for hazardous waste; 40 CFR268.40; Subpart D. D001: Ignitable characteristics of waste.

Table 1. Typical Data for ACCURUN HIV RNA Tri-Level Controls Kit.

ACCURUN HIV RNA Tri-Level Control	Analyte	Test	Typical Results
HIV RNA Tri-Level High Control	HIV RNA	Roche COBAS® AmpliPrep/COBAS® TaqMan® HIV RNA	271,487 copies/ml
HIV RNA Tri-Level Low Control	HIV RNA	Roche COBAS® AmpliPrep/COBAS® TaqMan® HIV RNA	344 copies/ml
HIV RNA Tri-Level Neg Control	HIV RNA	Roche COBAS® AmpliPrep/COBAS® TaqMan® HIV RNA	BLD

For assistance, contact SeraCare Technical Support at +1 508.244.6400.